

Claims Status

Claims 1, 2 and 21 are pending in this application. Claim 1 is amended herein. Support for the amendment to claim 1 can be found in the specification at pages 8-9.

Rejections Under 35 U.S.C. §112-Indefiniteness

Claims 1, 2 and 21 stand rejected as being indefinite for failing to recite steps of the claimed processes. In claim 1, the Examiner contends that the phrase “during a chromatography process” does not positively recite any chromatography process.

In claim 1 (part iii) and claim 24 (part iv), the Examiner contends that the step of “selecting fractions” is not described, and therefore renders the claim indefinite.

By this amendment, claim 1 has been amended for clarity to recite that the sample is purified using column chromatography. In our opinion, the Examiner's rejection based on the absence of the chromatography steps is not well-founded. The specification discloses that the column chromatography used to separate the prophenoloxidase *is not* limited to specific conditions. On page 8, the specification recites that the type and amount of chelating agent, the kind of resin, and the kind of solvent (*i.e.*, buffer) solution used for the chromatography can be easily determined by those skilled in the art.

In addition, the Examiner's attention is respectfully directed to the fact that one of the improvements of the claimed method over the prior art methods is that the fractions exhibiting phenoloxidase activity (*i.e.*, the fractions which contain *pro*phenoloxidase) can be separated using column chromatography including size-exclusion or ion exchange chromatography. Use of size-exclusion or ion exchange chromatography is superior to other techniques such as affinity chromatography, which is more expensive and labor intensive (see page 8, lines 19-22).

In the present invention, the composition that can detect β -1,3-glucan specifically down to 20 pg/ml can be purified by column chromatography without any further particular separation processes such as affinity chromatography.

Any size-exclusion/ion exchange column chromatography containing a compatible resin for separating prophenoloxidase is contemplated by the present invention. This is supported by the specification. For example, page 8-9 indicates that the resin is preferably a dextran- or vinyl-based resin such as Sephadex® or Toyopearl®. Both of these resins are appropriate for use in size-exclusion/ion exchange chromatography (see attached product information regarding Toyopearl® HW (used in Example 1, page 14, line 10; and product information regarding Sephadex® resins for ion-exchange chromatography-both at Tab A).

Moreover, as of the filing date, chromatography for separating proteins, including plasma proteins, was a well-known technique in the pertinent art. See for example, the attached abstract by Toribio et al. (Tab B).


Regarding the Examiner's rejection that the term "selecting fractions" is indefinite, the Examiner's attention is again directed to the specification at pages 10-11, which fully describes methods for determining phenoloxidase activity by β -1,3-glycan in the presence of calcium. It is respectfully asserted that where a technical process is well-known and readily practiced by those of ordinary skill in the art, the details of the process do not have to be recited in the claims. See *Shatterproof Glass Corp. v. Libbey Owens Ford Co.*, 758 F.2d 613 (Fed. Cir. 1985) ("The amount of detail required to be included in claims depends on the particular invention and the prior art, and is not to be viewed in the abstract but in conjunction with whether the specification is in compliance with the first paragraph of section 112"). See also *Hybritech Inc. v. Monoclonal Antibodies, Inc.*, 802 F.2d 1367 (Fed Cir 1986) (Holding that method claims to a "sandwich assay" were not indefinite since calculating the antibody affinities using this technique was known in the art).

Accordingly, it is submitted that the present claims are definite, as column chromatography was well known in the art as of the filing date. Withdrawal of this rejection is respectfully requested.

In view of the above amendment, applicant believes the pending application is in condition for allowance.

Dated: November 12, 2004

Respectfully submitted,

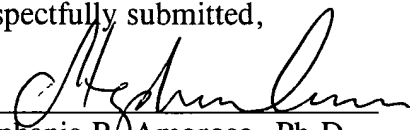By 
Stephanie R. Amoroso, Ph.D.

Registration No.: 51,401

DARBY & DARBY P.C.

P.O. Box 5257

New York, New York 10150-5257

(212) 527-7700

(212) 753-6237 (Fax)

Attorneys/Agents For Applicant



TOSOH BIOSCIENCE

TOSOH

Toyopearl® HW Media for Scale-up Size Exclusion Chromatography of Nucleic Acids and Proteins

CHEMISTRY OF TOYOPEARL RESINS

Toyopearl resins are hydrophilic macroporous packings for bioprocessing chromatography. These resins are semi-rigid, spherical beads synthesized by a copolymerization of ethylene glycol and methacrylate-type polymers. They are highly resistant to chemical attack and are not degraded by microbes. Toyopearl packings are available in the size ranges preferred for production chromatography and can be packed easily into columns up to the largest industrial size.

PROPERTIES OF TOYOPEARL HW MEDIA

Toyopearl HW media were developed for size exclusion (gel filtration) chromatography. There are five different products having increasing pore sizes and exclusion limits for proteins and water soluble polymers, as shown in **Table I**. **Table II** lists the features and benefits of Toyopearl HW resins.

Table I - Properties and Molecular Weight Separation Ranges for Toyopearl HW Resins


Toyopearl Resin	Particle Size (μm)	Pore Size (Å)	Molecular Weight of Sample		
			Globular Proteins	Dextrans	Polyethylene Glycols & Oxides
HW-40S HW-40F HW-40C	20 - 40 30 - 60 50 - 100	50	100 - 10,000	100 - 7,000	100 - 3,000
HW-50S HW-50F	20 - 40 30 - 60	125	500 - 80,000	500 - 20,000	100 - 18,000
HW-55S HW-55F	20 - 40 30 - 60	500	1,000 - 700,000	1,000 - 200,000	100 - 150,000
HW-65S HW-65F	20 - 40 30 - 60	1000	40,000 - 5,000,000	10,000 - 1,000,000	500 - 1,000,000
HW-75F	30 - 60	>1000	500,000 - 50,000,000	100,000 - 10,000,000	4,000 - 5,000,000

Table II - Features and Benefits of TSK-GEL HW Resins

- Constant packing volume over a wide range of salt concentrations
- High capacity for water soluble proteins
- High yield recovery using various aqueous eluents
- Stable from pH 2 to 12 for routine operation
Stable from pH 2 to 13 for cleaning
- Stable up to 3 bars of pressure
- Stable in eluents containing surfactants and organic solvents
- Can be autoclaved repeatedly in clean water at 120°C
- Bed stability
- Low cost chromatography
- Little loss or fouling
- Can be regenerated with base or acid
- Excellent flow rates, high linear velocities, can be scaled to very large columns
- Can be used for very hydrophobic proteins
- Easy to sanitize


FRACTIONATION RANGE

Figure 1 shows the fractionation range of the Toyopearl HW media. It can be seen that the product line covers the full range of normally encountered water soluble proteins. Also given is an indication of how Toyopearl HW resins compare with Sephadex® and Sephacryl® gels sold by Pharmacia.

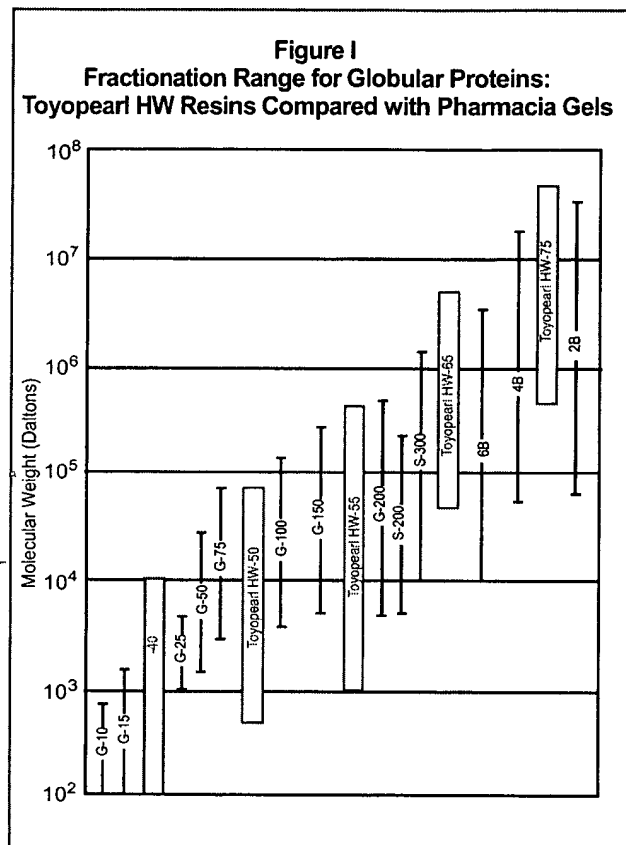


BED STABILITY IN VARIOUS SOLVENTS

Toyopearl Resins, due to their high degree of cross-linking, are quite resistant to changes in eluents as can be seen in **Figure 2**.

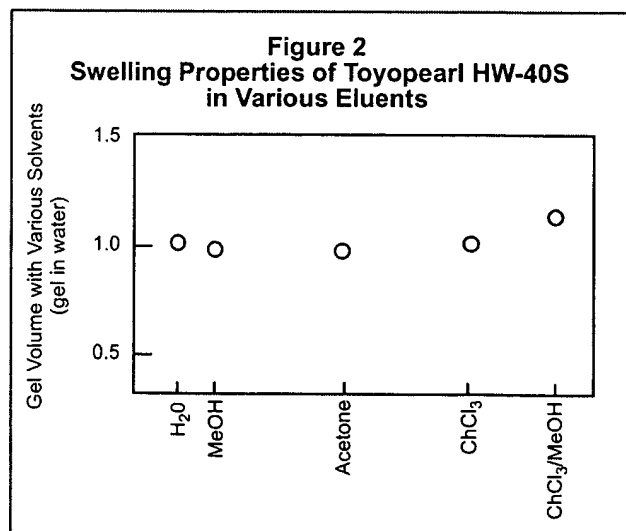


SIZE SORTING OF PROTEINS

The molecular weight of proteins can be estimated using size exclusion chromatography. **Figure 3** shows an example of this process using Toyopearl HW 55F in a chaotropic buffer.

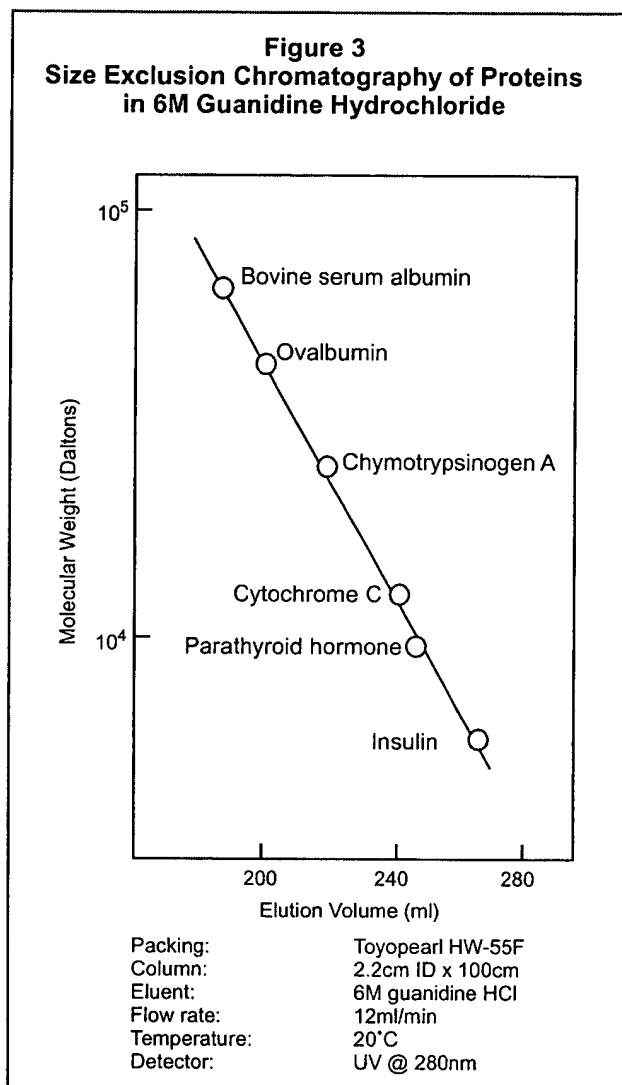


EFFECT OF SALT ON PROTEIN RETENTION

Protein retention on HW resins can change depending on the pI of the protein and the pH of the mobile phase. These effects can be minimized by the addition of a salt to the mobile phase as shown in **Figure 4**.

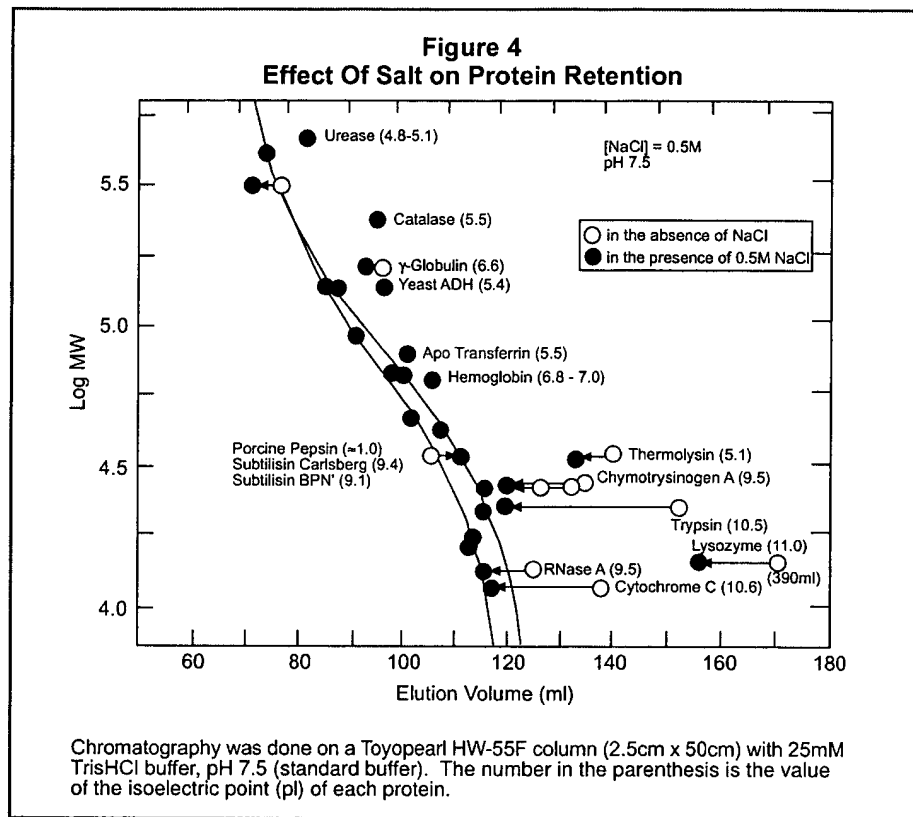


PRESSURE/FLOW CHARACTERISTICS

The superior rigidity of the Toyopearl Resin allows much higher flow rates than conventional gel products as illustrated in **Figure 5**.

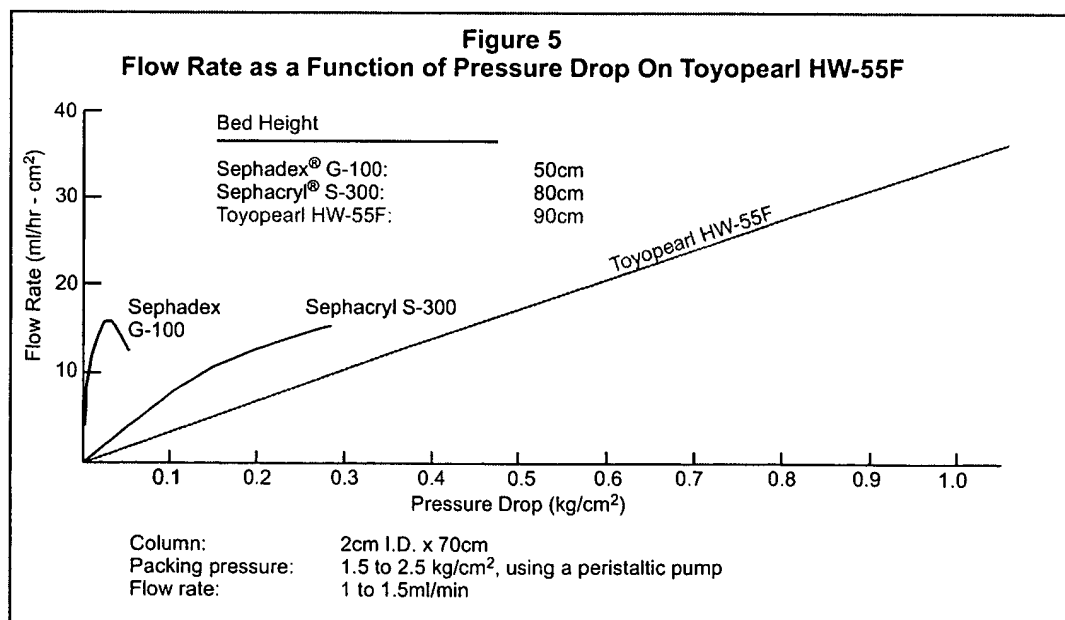

EFFECT OF PARTICLE SIZE ON RESOLUTION


In **Figure 6**, two different grades of Toyopearl HW-55 are compared --- an exercise often done in scale-up research and development. This figure shows that resolution is best with the S-grade (20 to 50 μ m), while the F-grade provides adequate performance. With these options, the most economical solution for scale-up chromatography can be chosen.

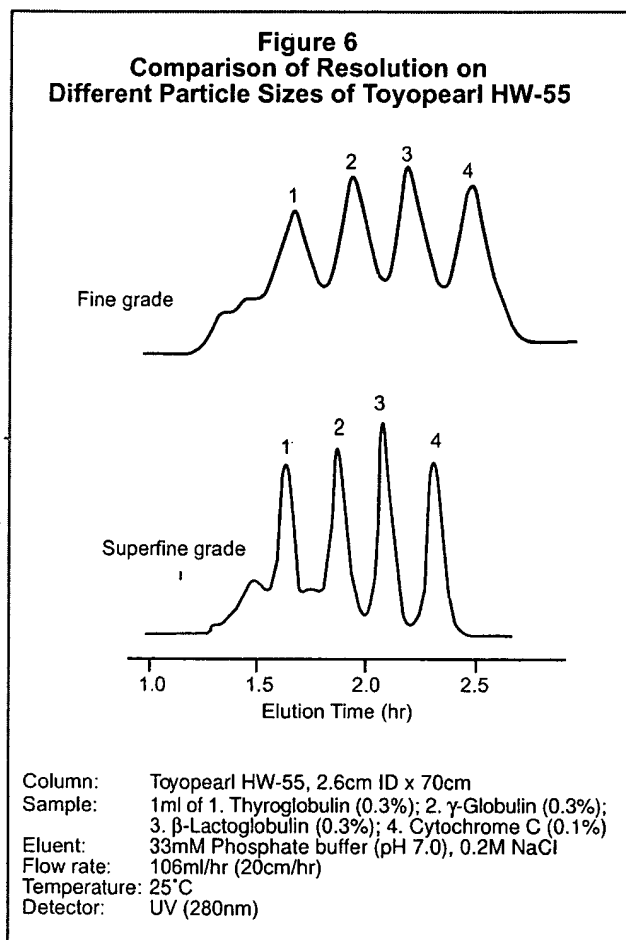

SCALING UP PROTEIN SEPARATIONS FOR HPLC

Due to surface and pore characteristics which are almost identical, Toyopearl HW resins behave very much as TSK-GEL HPLC columns as shown in **Figure 7**.


Figure 7
Separation of a Protein Mixture on
TSK-GEL G3000SW_{XL}, TSK-GEL G3000PW_{XL}
(5 μ m) HPLC Columns and Toyopearl HW-55F

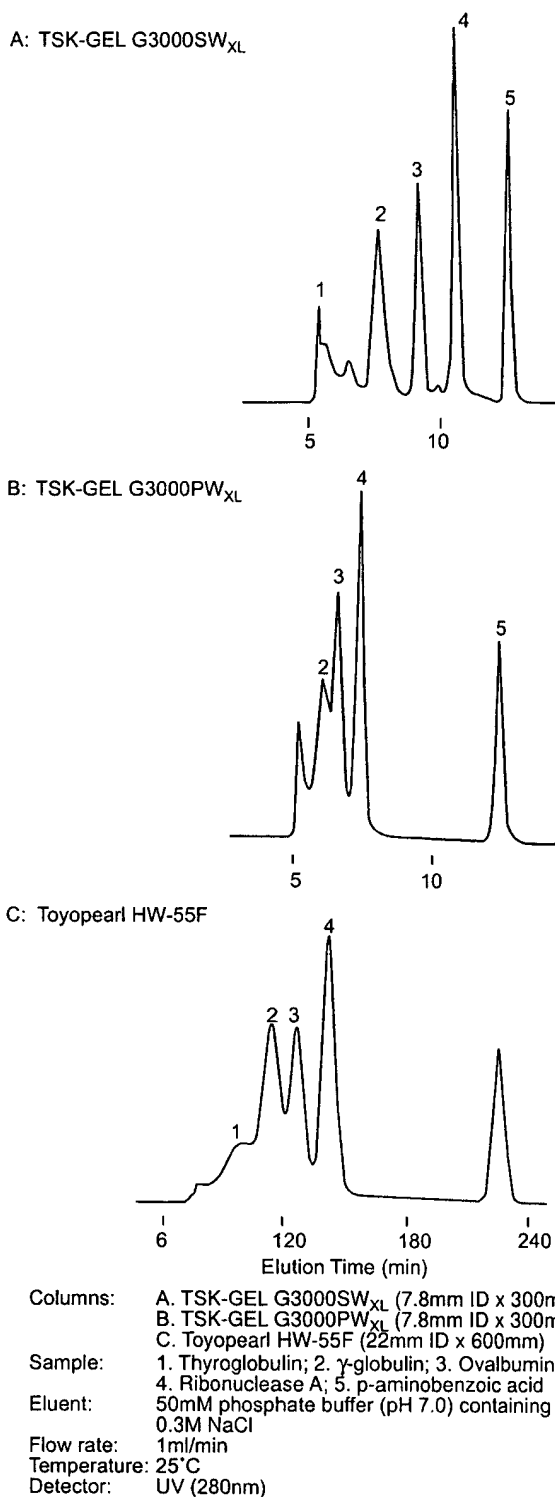

ADSORPTION CHROMATOGRAPHY EXAMPLES

SEPARATION OF NUCLEIC ACID-RELATED COMPOUNDS

Figure 8 shows a chromatogram of various nucleic acid related compounds obtained with Toyopearl HW-40S. This product is widely used for this kind of adsorption purification as well as for the separation of oligosaccharides.

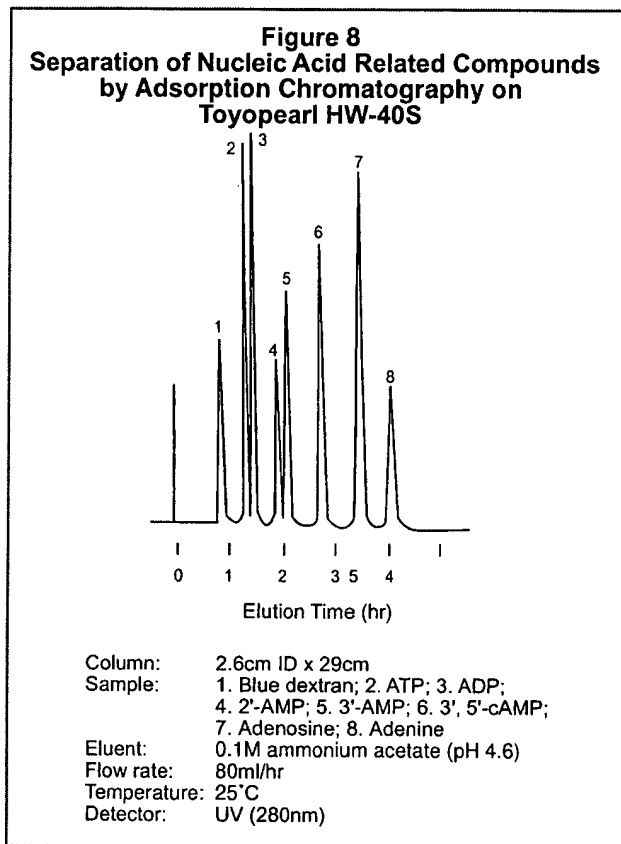

CHROMATOGRAPHY IN ORGANIC SOLVENTS

Figure 9 shows an example of purification of a plant-origin polyphenol in 50% methanol. Toyopearl HW-40 is being used in ways similar to ODS reverse phase packings or Sephadex LH-20 for the separation of polyphenolic substances.

ENZYME PURIFICATION EXAMPLES

PURIFICATION OF GLYCOPROTEINS

Membrane-bound enzymes, most of which are glycoproteins, tend to lose activity during conventional slow gel filtration. This process is due to excessive elution times as well as non-specific interaction between the glycoproteins and polysaccharide gel. High-speed elution with Toyopearl HW-55F solved this problem in the example shown in **Figure 10**. Recovery of enzymatic activity was up to 80%. Specific activity was increased as well, while chromatography with dextran gels destroyed 90% of the activity.

Figure 9
Elution Pattern of Polyphenolic Substances Obtained from Chestnut Galls


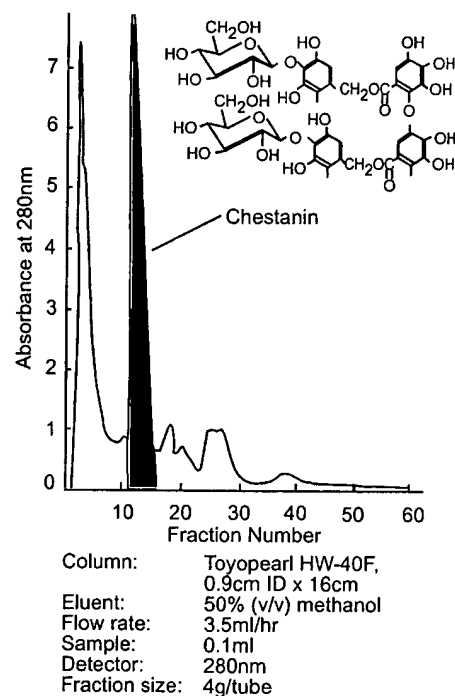


Figure 10
Chromatogram of Post-Proline Dipeptidyl Aminopeptidase from Pig Kidney

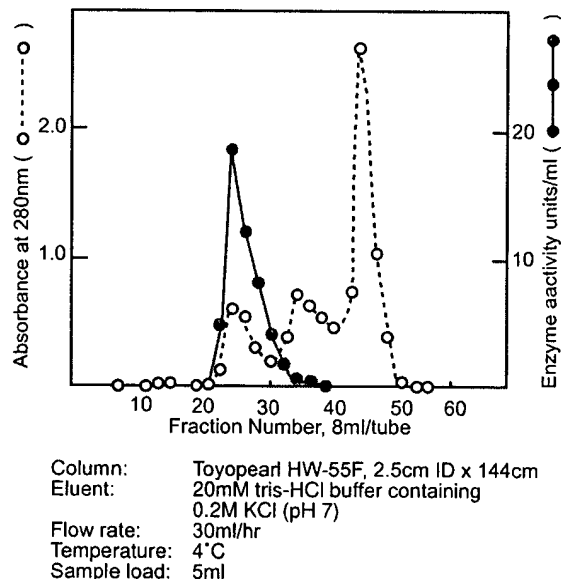

OTHER EXAMPLES OF ENZYME PURIFICATION

Figures 11 and 12 give two further examples of the use of Toyopearl HW resins for the recovery and purification of enzymes, one from a culture supernatant and the other from animal tissue.

Figure 11
Chromatogram of Extracellular Protease from *Thermus Calophilus* GK24



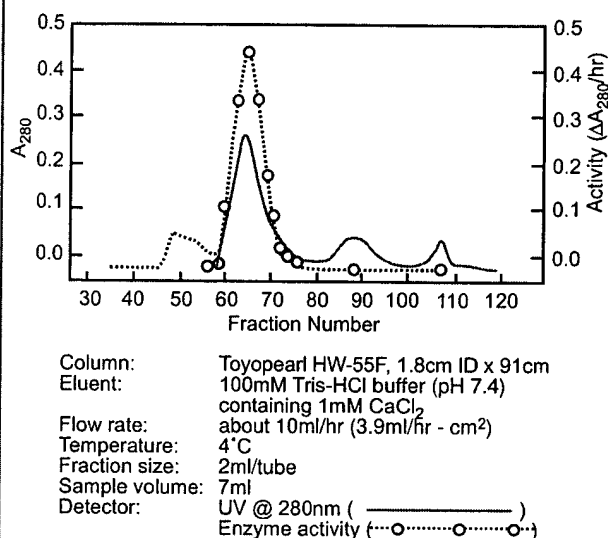


Figure 12
Purification of mRNA on Toyopearl HW-65F

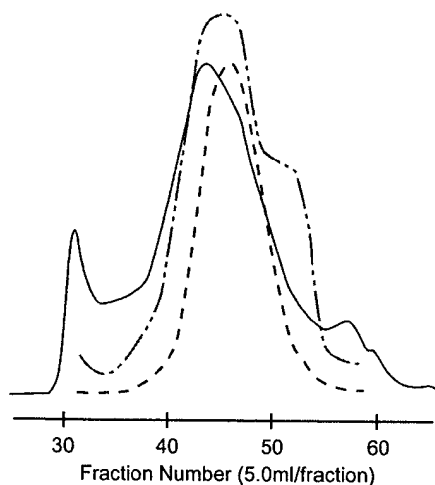


Toyopearl HW-50S resin has been used to help isolate the ubiquitin-histone conjugate uH2A from the unicellular ciliated protozoan *Tetrahymena pyriformis*. Figure 13 shows the separation of uH2A from the histone, H2A. The sole difference between these two components is a small polypeptide, ubiquitin (ca. 8,500Da).

Figure 13
Isolation of a Complex Protein Conjugate on Toyopearl HW-50S

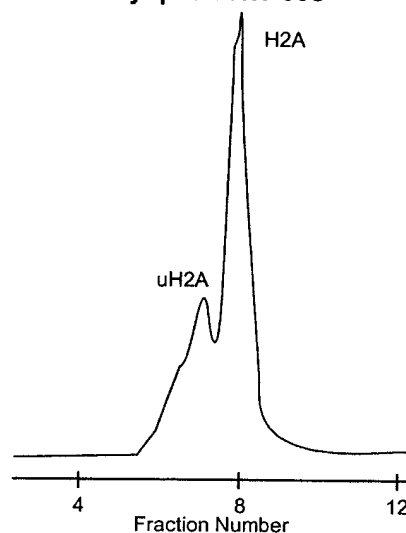

Gel filtration chromatography on a Toyopearl HW-75F column is also a sensitive and useful method for determining the purity of amylose specimens, and for demonstrating heterogeneity in molecular weight and branched structure. Figure 14 shows the elution profiles of amylose subfractions isolated from potato. Such substances as wheat and tapioca have also been purified by gel filtration.

Figure 14
Isolation of Potato Amylose Subfractions on Toyopearl HW-75F

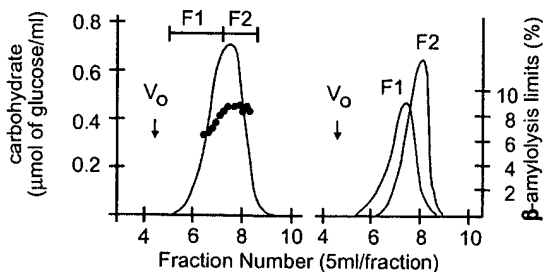

TOYOPEARL SEC RESINS:

Part #	Container size (ml)	Product description	Particle size (µm)	Exclusion limit (Da)
19809	150	HW-40S	20-40	3×10^3
07451	250			
14681	1000			
07967	5000			
19808	150	HW-40F	30-60	3×10^3
07448	500			
14682	1000			
07968	5000			
19807	150	HW-40C	50-100	3×10^3
07449	500			
14683	1000			
07969	5000			
19811	150	HW-50S	20-40	1.8×10^4
07455	250			
14684	1000			
08059	5000			
19810	150	HW-50F	30-60	1.8×10^4
07453	500			
14685	1000			
08060	5000			
19813	150	HW-55S	20-40	1.5×10^5
07459	250			
14686	1000			
08062	5000			
19812	150	HW-55F	30-60	1.5×10^5
07457	500			
14687	1000			
08063	5000			
19815	150	HW-65S	20-40	1×10^6
07467	250			
14688	1000			
08068	5000			
19814	150	HW-65F	30-60	1×10^6
07465	500			
14689	1000			
08069	5000			
19816	150	HW-75F	30-60	5×10^7
07469	500			
14691	1000			
08072	5000			

Conditions: Exclusion limits are +/- 30% and are determined using polyethylene glycol, polyethylene oxide, or dextran standards, as appropriate.

TOYOPEARL LABPAKs:

Part #	Container size (ml)	Product description	Particle size (µm)
19821	3 x 150	SECPAK LMW (HW-40F, HW-50F, HW-55F)	30-60
19819	3 x 150	SECPAK HMW (HW-55F, HW-65F, HW-75F)	30-60
19820	4 x 150	SECPAK HP (HW-40S, HW-50S, HW-55S, HW-65S)	20-40



TOSOH

TOSOH BIOSCIENCE

TOSOH BIOSCIENCE LLC
156 Keystone Drive
Montgomeryville, PA 18936-9637
Orders & Service: (800)366-4875
Fax: (215)283-5035
www.tosohbioscience.com
Orders & Service: sales&marketing@tosohbioscience.com


[Register for additional site benefits.](#) | [Login](#) | [Your Profile](#) | [Order Center](#) | [Search](#) | [Contact Us](#)

Aldrich

[Aldrich Handbook](#)[New Products](#)[Brand Info](#)[Aldrichimica Acta](#)[Aldrich Literature Request](#)[Equipment, Supplies & Books](#)[Tech Bulletins](#)

Technical Information Bulletin



Number AL-142

Ion-Exchange Resins and Related Polymeric Adsorbents

[Home](#) | [Amberlite® and Amberlyst® Resins](#) | [Dowex® Resins](#) | [Sephadex® Products](#) | [Merrifield Peptide Resins](#) | [Polyethylenimine-Based Products](#) | [Polyamide Resins](#)

III. SEPHADEX® PRODUCTS

These highly specialized gelfiltration and chromatographic media are composed of macroporous beads synthetically derived from the polysaccharide, dextran. The organic chains are covalently attached to the dextran backbone, giving a three-dimensional network having functional ionic groups attached by ether linkages to the glucose units of the polysaccharide chains.


Available forms include anion and cation exchangers, as well as gelfiltration resins, with degrees of porosity; bead sizes fall in discrete ranges between 20 and 300 μ . See Table XIII.

Table XII

Catalog Number	Description	Functionality	Bead size
Sephadex Ion-Exchange Resins			
27,124-1	Sephadex-CM C-25	carboxymethyl	40-60 μ
27,126-8	Sephadex-CM C-50	carboxymethyl	40-60 μ
27,127-6	Sephadex-DEAE A-25	2-(diethylamino)ethyl	40-60 μ
27,128-4	Sephadex-DEAE A-50	2-(diethylamino)ethyl	40-60 μ
27,129-2	Sephadex-QAE A-25	quaternary aminoethyl	40-60 μ
27,130-6	Sephadex-QAE A-50	quaternary aminoethyl	40-60 μ
27,131-4	Sephadex-SP C-25	sulfopropyl	40-60 μ
27,132-2	Sephadex-SP C-50	sulfopropyl	40-60 μ

Table XIII

Catalog Number	Description	Bead size
----------------	-------------	-----------

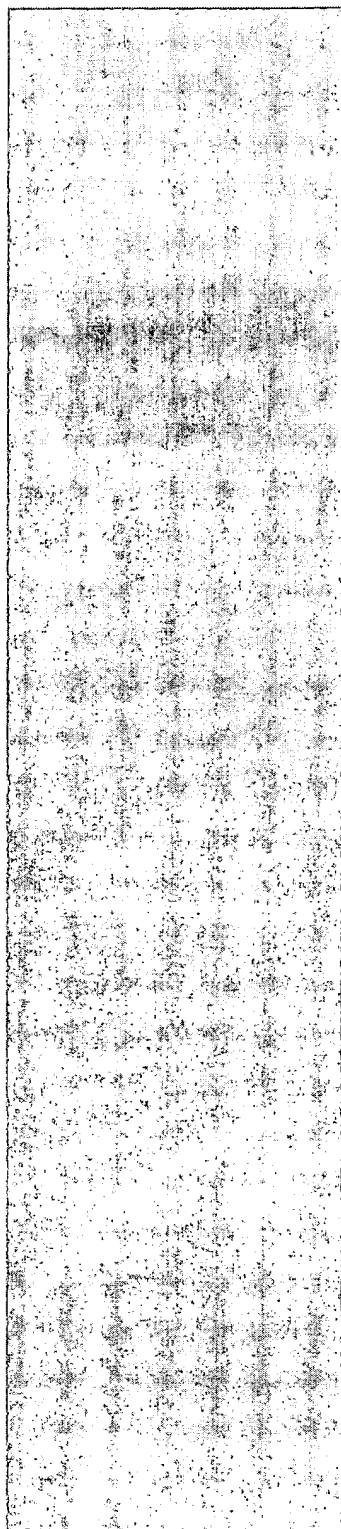




Sephadex Gel-Filtration Resins		
27,103-9	Sephadex G-10	40-120
27,104-7	Sephadex G-15	40-120
27,106-3	Sephadex G-25	20-50
27,107-1	Sephadex G-25	20-80
27,109-8	Sephadex G-25	50-150
27,110-1	Sephadex G-25	100-300
27,112-8	Sephadex G-50	20-50
27,113-6	Sephadex G-50	20-80
27,114-4	Sephadex G-50	50-150
27,115-2	Sephadex G-50	100-300
27,116-0	Sephadex G-75	20-50
27,117-9	Sephadex G-75	40-120
27,118-7	Sephadex G-100	20-50
27,119-5	Sephadex G-100	40-120
27,121-7	Sephadex G-150	40-120
27,123-3	Sephadex G-200	40-120

[Home](#) | [Amberlite® and Amberlyst® Resins](#) | [Dowex® Resins](#) | [Sephadex® Products](#) | [Merrifield Peptide Resins](#) | [Polyethylenimine-Based Products](#) | [Polyamide Resins](#)

Use of this web site constitutes your acceptance of the Site Use Terms

[help](#) | [privacy](#) | [technical library](#) | [search](#) | [home](#)
[terms and conditions](#) | [contract manufacturing](#)

© 2004 Sigma-Aldrich Co. Reproduction forbidden without permission.
Sigma-Aldrich brand products are sold exclusively through Sigma-Aldrich, Inc. Best viewed in IE5.

				National Library of Medicine 					
Entrez	PubMed	Nucleotide	Protein	Genome	Structure	OMIM	PMC	Journals	Books
Search <input type="text" value="PubMed"/>		for <input type="text"/>		<input type="button" value="Go"/>		<input type="button" value="Clear"/>			
Limits		Preview/Index		History		Clipboard		Details	
Display <input type="text" value="Abstract"/>		Show: <input type="text" value="20"/>		Sort <input type="text"/>		Send to <input type="text" value="Text"/>			
About Entrez		<input type="checkbox"/> 1: J Chromatogr B Biomed Appl. 1996 Sep 20;684(1-2):1-23. Related Article							
Text Version		Methods for chromatographic and electrophoretic separation and assay of NADP oxidoreductases.							
Entrez PubMed		Toribio F, Alhama J, Lopez-Barea J.							
Overview		Departamento de Bioquímica y Biología Molecular, Facultad de Veterinaria, Universidad de Córdoba, Spain.							
Help / FAQ									
Tutorial									
New/Noteworthy									
E-Utilities									
PubMed Services									
Journals Database									
MeSH Database									
Single Citation Matcher									
Batch Citation Matcher									
Clinical Queries									
LinkOut									
Cubby									
Related Resources									
Order Documents									
NLM Catalog									
NLM Gateway									
TOXNET									
Consumer Health									
Clinical Alerts									
ClinicalTrials.gov									
PubMed Central									
		<p>The different techniques described in purification protocols for pyridine nucleotide-dependent enzymes have been reviewed, covering mainly the papers published in the past six years. Chromatography was reported in 100% of reviewed papers and among the chromatographic techniques, affinity chromatography was the most used (ca. 92%), followed by ion-exchange chromatography (ca. 79%), size-exclusion chromatography (ca. 64%) and hydrophobic chromatography (ca. 24%). Other chromatographic techniques used infrequently. Each chromatographic technique has a different specific capacity and chemical selectivity and, therefore, the order of selection should be based on a precise knowledge of the nature of the sample and the amount of target enzyme that it contains. Analytical electrophoresis was used in about 9% of the reviewed papers, with denaturing polyacrylamide gel electrophoresis (PAGE) being the most widely used mode (ca. 92%), followed by native PAGE (ca. 48%). The use of isoelectric focusing was reported in 14% of the papers, preparative gel electrophoresis was used in only 8% of the cases. The use of electrophoretic techniques was reported in only a few papers. The use of continuous enzymatic activity assay methods (spectrophotometric) was found in most papers, while high-performance liquid chromatography-based methods (discontinuous assays) were reported in only 11% of the reviewed articles.</p>							
		<p>Publication Types:</p> <ul style="list-style-type: none">◦ Review◦ Review, Academic							
		PMID: 8906463 [PubMed - indexed for MEDLINE]							
Display <input type="text" value="Abstract"/>		Show: <input type="text" value="20"/>		Sort <input type="text"/>		Send to <input type="text" value="Text"/>			
Write to the Help Desk									